

NucleoSpin<sup>®</sup> RNA Virus
 NucleoSpin<sup>®</sup> RNA Virus F

MACHEREY-NAGEL

August 2022/Rev. 13



MACHEREY-NAGEL www.mn-net.com

## Contact MN

#### Germany and international

MACHEREY-NAGEL GmbH & Co. KG Valencienner Str. 11 · 52355 Düren · Germany Tel.: +49 24 21 969-0 Toll-free: 0800 26 16 000 (Germany only) E-mail: info@mn-net.com

#### **Technical Support Bioanalysis**

Tel.: +49 24 21 969-270 E-mail: tech-bio@mn-net.com

## USA

MACHEREY-NAGEL Inc. 924 Marcon Blvd. · Suite 102 · Allentown PA, 18109 · USA Toll-free: 888 321 6224 (MACH) E-mail: sales-us@mn-net.com

#### France

 MACHEREY-NAGEL SAS

 1, rue Gutenberg – BP135 · 67720 Hoerdt Cedex · France

 Tel.:
 +33 388 68 22 68

 E-mail:
 sales-fr@mn-net.com

 MACHEREY-NAGEL SAS (Société par Actions Simplifiée) au capital de 186600 €

 Siret 379 859 531 0020 · RCS Strasbourg B379859531 · N° intracommunautaire FR04 379 859 531


#### Switzerland

MACHEREY-NAGEL AG Hirsackerstr. 7 · 4702 Oensingen · Switzerland Tel.: +41 62 388 55 00 E-mail: sales-ch@mn-net.com

www.mn-net.com

# Viral RNA isolation

## Protocol at a glance (Rev. 13)

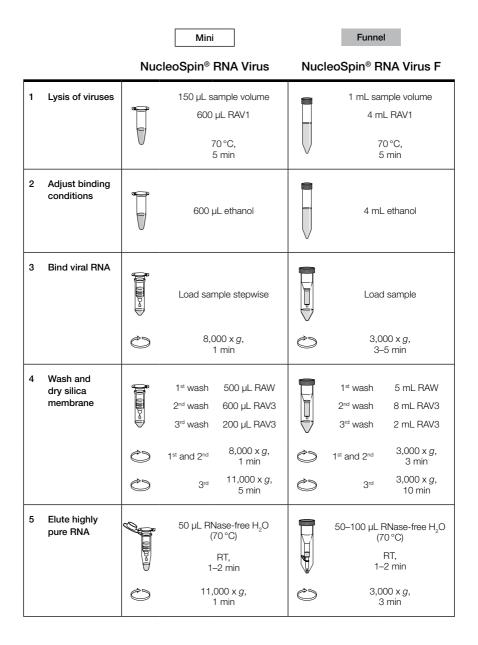



# Table of contents

1	1 Components	
	1.1 Kit contents	4
	1.2 Reagents, consumables, and equipment to be supplied by use	er 6
	1.3 About this user manual	6
2	2 Product description	7
	2.1 The basic principle	7
	2.2 Kit specifications	7
	2.3 Remarks regarding sample quality and preparation	8
	2.4 Remarks regarding elution	10
	2.5 Remarks regarding quality control	10
3	3 Storage conditions and preparation of working solutions	11
4	4 Safety instructions	13
	4.1 Disposal	13
5	5 NucleoSpin <sup>®</sup> RNA Virus protocols	14
	5.1 Viral RNA isolation from cell-free biological fluids	14
	5.2 Isolation of viral RNA and DNA from cell-free biological fluids	16
6	3 NucleoSpin <sup>®</sup> RNA Virus F protocols	17
	6.1 Viral RNA isolation from cell-free biological fluids	17
	6.2 Isolation of viral RNA and DNA from cell-free biological fluids	19
7	7 Appendix	20
	7.1 Troubleshooting	20
	7.2 Ordering information	21
	7.3 References	21
	7.4 Product use restriction / warranty	22

# 1 Components

## 1.1 Kit contents

	NucleoSpin <sup>®</sup> RNA Virus		
REF	10 preps 740956.10	50 preps 740956.50	250 preps 740956.250
Lysis Buffer RAV1	10 mL	35 mL	5 × 35 mL
Wash Buffer RAW	6 mL	30 mL	150 mL
Wash Buffer RAV3 (Concentrate)*	6 mL	12 mL	50 mL
RNase-free H <sub>2</sub> O	13 mL	13 mL	30 mL
Elution Buffer RE**	13 mL	13 mL	30 mL
Carrier RNA (lyophilized)	300 µg	1 mg	5 × 1 mg
NucleoSpin <sup>®</sup> RNA Virus Columns (dark blue rings, plus Collection Tubes)	10	50	250
Collection Tubes (2 mL)	30	150	750
User manual	1	1	1

<sup>\*</sup> For preparation of working solutions and storage conditions see section 3.

<sup>\*\*</sup> Composition of Elution Buffer RE: 5 mM Tris/HCl, pH 8.5

## Kit contents continued

	NucleoSpin <sup>®</sup> RNA Virus F
REF	25 preps 740958
Lysis Buffer RAV1	2 × 120 mL
Wash Buffer RAW	150 mL
Wash Buffer RAV3 (Concentrate)*	3 × 25 mL
RNase-free H <sub>2</sub> O	13 mL
Elution Buffer RE**	13 mL
Carrier RNA (lyophilized)	2 × 300 µg
NucleoSpin <sup>®</sup> RNA Virus F Columns (plus Collection Tubes)	25
Collection Tubes (50 mL)	25
Collection Tubes (0.5 mL)	25
User manual	1

<sup>\*</sup> For preparation of working solutions and storage conditions see section 3.

<sup>\*\*</sup> Composition of Elution Buffer RE: 5 mM Tris/HCl, pH 8.5

# 1.2 Reagents, consumables, and equipment to be supplied by user

Reagents

• 96-100 % ethanol

Consumables

- 1.5 mL microcentrifuge tubes (NucleoSpin<sup>®</sup> RNA Virus) or 50 mL tubes (NucleoSpin<sup>®</sup> RNA Virus F)
- Disposable tips

Equipment

- Manual pipettors
- Centrifuge for microcentrifuge tubes (NucleoSpin<sup>®</sup> RNA Virus) or 50 mL tubes (NucleoSpin<sup>®</sup> RNA Virus F)
- Vortex mixer
- Heating block or water bath for 70 °C incubation
- Personal protection equipment (lab coat, gloves, goggles)

## 1.3 About this user manual

It is strongly recommended reading the detailed protocol sections of this user manual if the **NucleoSpin<sup>®</sup> RNA Virus / RNA Virus F** kit is used for the first time. Experienced users, however, may refer to the Protocol at a glance instead. The Protocol at a glance is designed to be used only as a supplemental tool for quick referencing while performing the purification procedure.

All technical literature is available on the internet at *www.mn-net.com*.

Please contact Technical Service regarding information about changes of the current user manual compared to previous revisions.

## 2 Product description

## 2.1 The basic principle

With the **NucleoSpin<sup>®</sup> RNA Virus** method, RNA viruses are lysed quickly and efficiently by Lysis Buffer RAV1 which is a highly concentrated solution of GITC. DNA viruses (e.g., HBV) are usually more difficult to lyse and require Proteinase K digestion (see support protocol, section 5.2). Lysis buffer and ethanol create appropriate conditions for binding of nucleic acids to the silica membrane of the **NucleoSpin<sup>®</sup> RNA Virus Columns**. Carrier RNA improves binding and recovery of low-concentrated viral RNA. Contaminations (potential PCR inhibitors) like salts, metabolites and soluble macromolecular cellular components are removed in simple washing steps with ethanolic buffers **RAW** and **RAV3**. The nucleic acids can be eluted in low salt buffer or water and are ready-for-use in subsequent reactions.

## 2.2 Kit specifications

**NucleoSpin<sup>®</sup> RNA Virus / Virus F** kits are designed for the rapid preparation of highly pure viral nucleic acids (e.g., HCV, HIV, CMV) from fluid biological samples, for example plasma, serum, urine, swab wash solution, liquefied sputum, but not undiluted whole blood (see remarks in section 2.1).

- No cross contamination due to closed systems.
- The NucleoSpin<sup>®</sup> RNA Virus kit works with 150 µL serum, NucleoSpin<sup>®</sup> RNA Virus F funnel columns allow the processing of 1 mL serum.
- The funnel column of the NucleoSpin<sup>®</sup> RNA Virus F kit allows a high loading capacity as well as a simultaneously small elution volume.
- The prepared nucleic acids are suitable for applications like automated fluorescent DNA sequencing, RT-PCR\*, or any kind of enzymatic reaction. The detection limit for certain viruses depends on the individual procedures, for example in-house nested (RT-) PCR. We highly recommend using internal (low-copy) standards as well as positive and negative controls to monitor the purification, amplification, and detection processes.
- Carrier RNA (poly(-A) RNA: poly(A) potassium salt, prepared from ADP with polynucleotide phosphorylase) is included for optimal performance. Carrier RNA enhances binding of viral nucleic acids to the silica membrane and reduces the risk of viral RNA degradation. Please note that eluates of the NucleoSpin® RNA Virus kit contain both viral nucleic acids and Carrier RNA with amounts of Carrier RNA that may exceed the amount of viral nucleic acids. Therefore, it is not possible to quantify the nucleic acids isolated with the kit by photometric or fluorometric methods when using the carrier. Thus, other methods for quantification such as specific quantitative PCR or RT-PCR systems are recommended. Furthermore, Carrier RNA may inhibit PCR reactions. The amount of added Carrier RNA may thus be carefully optimized depending on the individual PCR system used.

Table 1: Kit specifications at a glance			
Parameter	NucleoSpin <sup>®</sup> RNA Virus	NucleoSpin <sup>®</sup> RNA Virus F	
Technology	Silica membrane technology	Silica membrane technology	
Format	Mini spin columns	Funnel columns	
Sample material	$\leq$ 150 µL serum, plasma, cell- free biological fluids	≤ 1 mL serum, plasma, cell-free biological fluids	
Fragment size	100 b–approx. 50 kb	100 b-approx. 50 kb	
Typical recovery rates	> 90 %	> 90 %	
Typical analysis limit	30–60 cp/mL*	30–60 cp/mL*	
Elution volume	50 μL	50–100 μL	
Preparation time	30 min/4–6 preps	45 min/2-4 preps	
Binding capacity	40 µg	30 µg	

## 2.3 Remarks regarding sample quality and preparation

### Liquid samples

Biological fluids or semi-fluid samples can be processed (e.g., serum, urine, or bronchoalveolar lavage). For successful nucleic acid purification it is important to obtain a homogeneous, clear, and non-viscous sample before loading onto the NucleoSpin<sup>®</sup> RNA Virus Column. Therefore, check all samples (especially old or frozen ones) for presence of precipitates. Precipitates remaining after lysis with Buffer RAV1 can be removed by centrifugation. Avoid clearing samples before lysis, because viruses of interest may be associated with particles or aggregates. Incubation with Buffer RAV1 can be prolonged in order to dissolve and digest residual cell structures, precipitates and virus particles. RNA, however, is sensitive and prolonged incubation may cause decreased yields.

#### Solid samples (tissue samples, stool samples)

Prepare a 10 % (w/v) suspension of tissue in buffer (e.g., PBS) using commercial homogenization tools (rotor-stator or bead-based homogenization tools, etc.). Centrifuge the suspension in order to remove particles. Use the clear particle-free supernatant for further processing.

<sup>\*</sup> Nested PCR

#### Swab material

Incubate swabs with moderate shaking for 30 min in a suitable buffer (e.g., PBS or liquid medium) to release sample material from the swab. Swab heads should be completely submerged in buffer (usually requires  $400-500 \mu$ L buffer). Transfer 150  $\mu$ L of the swab wash solution to a suitable container (e.g. reaction tube) and proceed with the sample lysis according to the standard protocol (section 5).

#### Viscous biological fluids (sputum, bronchoalveolar lavage, etc.)

Viscous biological samples, such as sputum and bronchoalveolar lavage, should be liquefied before starting the extraction protocol. Transfer 150  $\mu$ L sample to a suitable reaction tube. Add 600  $\mu$ L Buffer RAV1 containing Carrier RNA to the sample and incubate at 70 °C for 10 min with moderate shaking. Allow the sample to cool down, check if the sample is liquefied and proceed with adjusting the binding conditions according to the standard protocol (section 5).

In case the sample is not liquefied after the heat incubation restart the procedure following recommended guidelines (e.g., CDC guidelines for processing sputum samples:

https://www.cdc.gov/coronavirus/2019-ncov/downloads/processing-sputum-specimens.pdf). Transfer 150 µL of the liquefied sample to a suitable reaction tube and proceed with the standard protocol (section 5) starting with the sample lysis.

#### Blood samples

Processing of blood samples is possible if using blood diluted with PBS buffer. Using undiluted blood may cause clogging of the silica membrane of the NucleoSpin<sup>®</sup> Virus Binding Plate. The amount of PBS buffer added to blood samples has to be optimized for the individual organism. As a rule of thumb we recommend to start with 50  $\mu$ L blood diluted with 50  $\mu$ L PBS buffer.

#### Proteinase K treatment

Addition of Proteinase K solution (see ordering information) is necessary for the isolation of viral DNA or simultaneous viral RNA/DNA isolation. For isolation of viral RNA, Proteinase K treatment is usually not required. Proteinase K treatment is recommended for viral RNA isolation when viscous samples have to be processed (e.g., sputum samples).

#### Sample lysis

For isolation of viral RNA in general a lysis of samples in Buffer RAV1 for 10 min at room temperature will be sufficient. For isolation of viral RNA from viscous samples, for example sputum or supernatants of tissue suspensions or stool samples, a lysis at 70 °C may be required. For simultaneous isolation of viral RNA and DNA, incubation time (e.g., 5–15 min), and temperature (e.g., RT, 56 °C, or 70 °C) should be optimized and adjusted to the sample material used.

## 2.4 Remarks regarding elution

- Pure nucleic acids are finally eluted under low ionic strength conditions with RNase-free H<sub>2</sub>O (pH about 7 – 8) or slightly alkaline Buffer RE (5 mM Tris-HCl, pH 8.5).
- Elution can be performed in a single step with water/elution buffer as indicated in the
  protocol, obtaining at least 80% of the bound nucleic acids. To improve sensitivity,
  this eluate can be used in a second elution step increasing the efficiency of elution and
  concentration of viral nucleic acids slightly. Alternatively, a second elution step can be
  performed with an additional volume of water/ elution buffer releasing practically all
  bound nucleic acids but resulting in a lower concentrated, combined eluate.
- RNA should be eluted with the water supplied and DNA with Elution Buffer RE. Buffer RE provides better storage conditions for DNA. To elute both types of nucleic acids together, use the pH proofed (pH 6 8), RNase-free H<sub>2</sub>O preheated to 70 °C.

## 2.5 Remarks regarding quality control

 Buffers and NucleoSpin<sup>®</sup> RNA Virus / RNA Virus F Columns have been tested with rRNA and MS2 phage RNA. The absence of RNases, and the yield and efficiency of purification have been investigated with RT-PCR.

# 3 Storage conditions and preparation of working solutions

Attention: Buffers RAV1 and RAW contain chaotropic salts. Wear gloves and goggles!

CAUTION: Buffer RAV1 contains guanidinium thiocyanate and Buffer RAW contains guanidine hydrochloride which can form highly reactive compounds when combined with bleach (sodium hypochlorite). DO NOT add bleach or acidic solutions directly to the sample-preparation waste.

- All kit components can be stored at room temperature (15-25 °C) and are stable until: see package label.
- Carrier RNA has a limited shelf life in Buffer RAV1. For this reason, some kits contain several bottles of lyophilized Carrier RNA that should be used successively as required, to avoid degradation of Carrier RNA.

<u>Note:</u> Due to the production procedure and the small amount of Carrier RNA contained in the vial, the carrier may hardly be visible in the vial.

 Before use, add 1 mL Lysis Buffer RAV1 to the Carrier RNA tube. Dissolve the RNA and transfer it back to the RAV1 bottle.

Storage of Carrier RNA in Buffer RAV1:

Lysis Buffer RAV1 including Carrier RNA can be stored at room temperature for 1 – 2 weeks. Storage at room temperature prevents salt precipitation.
 Lysis Buffer RAV1 including Carrier RNA can be stored at 4 °C for up to 4 weeks or aliquoted and stored at -20 °C for longer periods. Carrier-RNA dissolved in RAV1 and stored at -20 °C is stable for at least one year. Storage at 4 °C or below may cause salt precipitation. Therefore, the mixture must be preheated at 40–60 °C for a maximum of 5 min in order to redissolve salts.
 Do not warm Buffer RAV1 containing Carrier RNA more than 4 times! Frequent

warming, temperatures > 80 °C, and extended heat incubation will accelerate the degradation of Carrier RNA. This leads to reduced recovery of viral RNA and eventually false negative RT-PCR results, in particular if low-titer samples are used.

Before starting any NucleoSpin<sup>®</sup> RNA Virus/RNA Virus F protocol, prepare the following:

• Wash Buffer RAV3: Add the indicated volume of ethanol (96–100%) to Wash Buffer RAV3 Concentrate. Mark the label of the bottle to indicate that the ethanol is added. Store Wash Buffer RAV3 at room temperature for up to one year.

	NucleoSpin <sup>®</sup> RNA Virus		
REF	10 preps	50 preps	250 preps
	740956.10	740956.50	740956.250
Wash Buffer RAV3	6 mL	12 mL	50 mL
(Concentrate)	Add 24 mL ethanol	Add 48 mL ethanol	Add 200 mL ethanol
	NucleoSpin <sup>®</sup> RNA Virus F		
REF	25 preps 740958		
Wash Buffer RAV3	3 × 25 mL		
(Concentrate)	Add 100 mL ethanol to each vial		

# 4 Safety instructions

When working with the **NucleoSpin<sup>®</sup> RNA Virus** kit wear suitable protective clothing (e.g., lab coat, disposable gloves, and protective goggles). For more information consult the appropriate Material Safety Data Sheets (MSDS available online at *http://www.mn-net.com/msds*).



Caution: Guanidine hydrochloride in Buffer RAW guanidinium thiocyanate in Buffer RAV1 can form highly reactive compounds when combined with bleach! Thus, do not add bleach or acidic solutions directly to the sample preparation waste.

The waste generated with the **NucleoSpin® RNA** Virus kit has not been tested for residual infectious material. A contamination of the liquid waste with residual infectious material is highly unlikely due to strong denaturing lysis buffer treatment but it cannot be excluded completely. Therefore, liquid waste must be considered infectious and should be handled and discarded according local safety regulations.

## 4.1 Disposal


Dispose hazardous, infectious or biologically contaminated materials in a safe and acceptable manner and in accordance with all local and regulatory requirements.


# 5 NucleoSpin<sup>®</sup> RNA Virus protocols

## 5.1 Viral RNA isolation from cell-free biological fluids

Before starting the preparation:

- Check if Lysis Buffer RAV1 and Wash Buffer RAV3 were prepared according to section 3.
- Preheat an aliquot of Elution Buffer RE / RNase-free H<sub>2</sub>O to 70 °C.
- Make sure to prepare the sample according to the remarks for sample preparation (section 2.3).

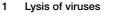

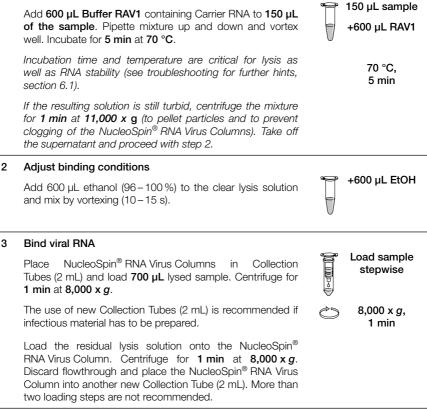

## 4 Wash and dry silica membrane

## 1 st wash

Add **500 μL Buffer RAW** to the NucleoSpin<sup>®</sup> RNA Virus Column. Centrifuge for **1 min** at **8,000 x g**. Discard flowthrough.

This washing step removes contaminants and PCR inhibitors.

## 2<sup>nd</sup> wash

Add **600 \muL Buffer RAV3** to the NucleoSpin<sup>®</sup> RNA Virus Column. Centrifuge for **1 min** at **8,000 x** *g*. Discard flowthrough with Collection Tube.

## 3rd wash

Place the NucleoSpin<sup>®</sup> RNA Virus Column in a new Collection Tube (2 mL) and add **200 µL Buffer RAV3**. Centrifuge for **2–5 min** at **11,000 x** *g* to remove ethanolic Buffer RAV3 completely.

<u>Optional:</u> Residual Buffer RAV3 may inhibit subsequent reactions. Therefore, for subsequent reactions which are extremely ethanol-sensitive, we recommend repeating the centrifugation with a new Collection Tube (2 mL). Or alternatively, incubate the NucleoSpin<sup>®</sup> RNA Virus Columns for 1 min at 70 °C to remove any remaining traces of ethanol.

#### 5 Elute viral RNA

Place the NucleoSpin<sup>®</sup> RNA Virus Column into a new, sterile 1.5 mL microcentrifuge tube (not provided). Add 50  $\mu$ L RNase-free H<sub>2</sub>O (preheated to 70 °C) and incubate for 1–2 min. Centrifuge for 1 min at 11,000 x g.

To elute viral DNA which was prepared according to the support protocol 5.2, we recommend using Buffer RE, preheated to 70  $^{\circ}$ C (also see section 2.4).

+600 μL RAV3

+500 µL RAW

8,000 x q,

1 min

8,000 x *g*, 1 min

+200 µL RAV3

11,000 x g, 5 min

# 5.2 Isolation of viral RNA and DNA from cell-free biological fluids

This protocol is recommended for the purification of viral RNA and viral DNA for all types of DNA viruses like HBV and CMV from small samples of up to 150  $\mu$ L.

## Before starting the preparation:

- Check if Lysis Buffer RAV1 and Wash Buffer RAV3 were prepared according to section 3.
- Check if Proteinase K solution (not included, see ordering information) was prepared.
- Preheat an aliquot of Elution Buffer RE / RNase-free H<sub>2</sub>O to 70 °C.
- Make sure to prepare the sample according to the remarks for sample preparation (section 2.3).

#### 1 Lysis of viruses

Add 600  $\mu$ L Buffer RAV1 containing Carrier RNA to 150  $\mu$ L of the sample. Add 20  $\mu$ L Proteinase K (20 mg/mL stock solution), to the lysis mixture. Pipette mixture up and down and vortex for 10–15 s. Incubate for 5 min at 70 °C.

Incubation time and temperature are critical for lysis as well as RNA stability (see troubleshooting).

Proteinase K is not included in this kit, but can be ordered separately (see ordering information). For the isolation of viral DNA and genomic DNA from other matrices (**not cell-free**) like blood we recommend the **NucleoSpin<sup>®</sup> Blood** or **NucleoSpin<sup>®</sup> Tissue** kit (see ordering information).

If the resulting solution is still turbid, centrifuge the mixture for **1 min** at **11,000 x g** to pellet particles and to prevent clogging of the NucleoSpin<sup>®</sup> RNA Virus Columns. Take off the supernatant and continue with step 2 of protocol 5.1. 150 μL sample +600 μL RAV1 +20 μL Prot. K

> 70 °C, 5 min

# 6 NucleoSpin<sup>®</sup> RNA Virus F protocols

## 6.1 Viral RNA isolation from cell-free biological fluids

## Before starting the preparation:

- Check if Lysis Buffer RAV1 and Wash Buffer RAV3 were prepared according to section 3.
- Preheat an aliquot of Elution Buffer RE/RNase-free H<sub>2</sub>O to 70 °C.
- Make sure to prepare the sample according to the remarks for sample preparation (section 2.3).

#### 1 Lysis of viruses

Add **4 mL Buffer RAV1** containing Carrier RNA to **1 mL of the sample**. Pipette mixture up and down and vortex well. Incubate for **5 min** at **70 °C**.

Incubation time and temperature are critical for lysis as well as RNA stability (see troubleshooting for further hints).

If the resulting solution is still turbid, centrifuge the mixture for **1 min** at **11,000 x g** (to pellet particles and to prevent clogging of the NucleoSpin<sup>®</sup> RNA Virus F Columns). Take off the supernatant and proceed with step 2. 1 mL sample +4 mL RAV1

> 70 °C, 5 min

+4 mL EtOH

## 2 Adjust binding conditions

Add **4 mL ethanol** (96-100%) to the clear lysis solution and mix by vortexing (10-15 s).

## 3 Bind viral RNA

Take the NucleoSpin<sup>®</sup> RNA Virus F Column placed in a Collection Tube and load lysed sample. Centrifuge for  $3-5 \min$  at  $3,000 \times g$ .

The use of new 50 mL tubes for every step is recommended if infectious material has to be prepared. This avoids crosscontamination and contamination of centrifuge units. For non-infectious samples, we recommend discarding the flowthrough and reusing the 50 mL tube for loading and washing steps. Additional 50 mL tubes have to be ordered separately.

The maximum loading capacity of the NucleoSpin<sup>®</sup> RNA Virus F Column is about 10 mL in order to work crosscontamination free. If more sample has to be loaded discard flowthrough and put the NucleoSpin<sup>®</sup> RNA Virus F Column into a new 50 mL tube. Load sample

3,000 x g, 3–5 min Load the residual lysis solution onto the NucleoSpin<sup>®</sup> RNA Virus F Column. Centrifuge for 3-5 min at  $3,000 \times \text{g}$ . Discard flowthrough and place the NucleoSpin<sup>®</sup> RNA Virus F Column into another new 50 mL tube. More than two loading steps are not recommended.

4 Wash and dry silica membrane

## 1 st wash

Add **5 mL Buffer RAW** to the NucleoSpin<sup>®</sup> RNA Virus F Column. Centrifuge for **3 min** at **3,000 x** *g*. Discard flowthrough.

This washing step removes contaminants and PCR inhibitors.

#### 2<sup>nd</sup> wash

Add **8 mL Buffer RAV3** to the NucleoSpin<sup>®</sup> RNA Virus F Column. Centrifuge for **3 min** at **3,000 x** *g*. Discard flowthrough with Collection Tube.


#### 3<sup>rd</sup> wash

Place the NucleoSpin<sup>®</sup> RNA Virus F Column in a new Collection Tube (50 mL) and add **2 mL Buffer RAV3**. Centrifuge for **10 min** at **3,000 x** *g* to remove ethanolic Buffer RAV3 completely.

<u>Optional:</u> Residual Buffer RAV3 may inhibit subsequent reactions. Therefore, for subsequent reactions, which are extremely ethanol-sensitive, we recommend repeating the centrifugation with a new Collection Tube (50 mL). Or alternatively, incubate the NucleoSpin<sup>®</sup> Virus F Columns for 1 min at 70 °C to remove any remaining traces of ethanol.

## 5 Elute viral RNA

Attach the supplied Collection Tube (0.5 mL) with the adaptor to the NucleoSpin<sup>®</sup> RNA Virus F Column. Place the assembly in a 50 mL tube (not provided). Add  $50-100 \,\mu$ L RNase-free H<sub>2</sub>O (preheated to 70 °C) and incubate for  $1-2 \,min$  at room temperature. Centrifuge for  $3 \,min$  at 3,000 x g.

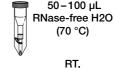

1 – 2 min

3,000 x g, 3 min



+5 mL RAW

3,000 x g,

3 min

+2 mL RAV3 3.000 x *a*.

3,000 x g 10 min

# 6.2 Isolation of viral RNA and DNA from cell-free biological fluids

This protocol is recommended for the purification of viral RNA and viral DNA for all types of DNA viruses like HBV and CMV for samples of up to 1 mL.

## Before starting the preparation:

- Check if Lysis Buffer RAV1 and Wash Buffer RAV3 were prepared according to section 3.
- Check if Proteinase K solution (not included, see ordering information) was prepared.
- Preheat an aliquot of Elution Buffer RE/RNase-free H<sub>2</sub>O to 70 °C.
- Make sure to prepare the sample according to the remarks for sample preparation (section 2.3).

#### 1 Lysis of viruses

Add 4 mL Buffer RAV1 containing Carrier RNA to 1 mL of the fluid sample. Add 133  $\mu$ L Proteinase K (20 mg/mL stock solution), to the lysis mixture. Pipette mixture up and down and vortex for 10–15 s. Incubate for 5 min at 70 °C.

Incubation time and temperature are critical for lysis as well as RNA stability (see troubleshooting).

Proteinase K is not included in this kit, but can be ordered separately (see ordering information). For the isolation of viral DNA and genomic DNA from other matrices (not cell-free) like blood we recommend the **NucleoSpin<sup>®</sup> Blood** or **NucleoSpin<sup>®</sup> Tissue** kit (see ordering information).

If the resulting solution is still turbid, centrifuge the mixture for **1 min** at **11,000 x g** to pellet particles and to prevent clogging of the NucleoSpin<sup>®</sup> RNA Virus Columns. Take off the supernatant and continue with step 2 of protocol 5.1. 1 mL sample
 +4 mL RAV1
 +133 μL Prot. K

70 °C, 5 min

# 7 Appendix

## 7.1 Troubleshooting

Problem	Possible cause and suggestions		
	Problems with Carrier RNA		
	Carrier RNA not added.		
	<ul> <li>See remarks concerning storage of Buffer RAV1 with Carrier RNA (section 3).</li> </ul>		
Small	Proteinase K digestion may be necessary		
amounts or no viral	<ul> <li>Use and compare protocols with and without Proteinase K digestion or prolong incubation time to 10 min.</li> </ul>		
nucleic acids in the eluate	Viral nucleic acids degraded		
	<ul> <li>Samples should be processed immediately. If necessary, add RNase inhibitor to the sample and ensure appropriate storage conditions up to the processing.</li> </ul>		
	<ul> <li>Check that all buffers have been prepared and stored correctly. If in doubt, use new aliquots of Buffer RAV1, Carrier RNA and Elution Buffer RE.</li> </ul>		
	Reduced sensitivity		
	• Change the volume of eluate added to the PCR/RT-PCR.		
Problems with subsequent detection	<ul> <li>Incubation time and temperature are critical for lysis as well as RNA stability. For sensitive RNA preparations, incubation at room temperature is sufficient without significant loss of sensitivity. For parallel isolation of viral RNA and DNA incubation time (5–15 min) and temperature (RT/56 °C/72 °C) may be adapted in order to get optimal recovery rates for both species.</li> </ul>		
	Ethanol carry-over		
	<ul> <li>Prolong centrifugation steps in order to remove Buffer RAV3 completely.</li> </ul>		
General	Clogged membrane		
problems	<ul> <li>Centrifuge plasma lysate before the addition of ethanol and subsequent loading onto the corresponding NucleoSpin<sup>®</sup> RNA Virus Columns.</li> </ul>		

## 7.2 Ordering information

Product	REF	Pack of
NucleoSpin <sup>®</sup> RNA Virus	740956.10/.50/.250	10/50/250
NucleoSpin <sup>®</sup> RNA Virus F	740958	25
NucleoSpin <sup>®</sup> Funnel Columns	740959	30 sets
Proteinase K	740506	100 mg
rDNase Set (Recombinant DNase and Reaction Buffer for rDNase; sufficient for 50 mini preps)	740963	1 set
NucleoSpin <sup>®</sup> Dx Virus	740895.50/.250	50/250
NucleoSpin <sup>®</sup> Blood	740951.10/.50/.250	10/50/250
NucleoSpin <sup>®</sup> RNA Blood	740200.10/.50	10/50
NucleoSpin <sup>®</sup> Tissue	740952.10/.50/.250	10/50/250
Collection Tubes (2 mL)	740600	1000

Visit www.mn-net.com for more detailed product information.

## 7.3 References

M.L. Villahermosa, M.Thomson, E.Vazques de Parga, M.T. Cuevas, G. Contreas, L.Perez-Alvarez, E.Delgado, N.Manjon, L.Medrano and R.Najera. Improved Conditions for Extraction and Amplification of Human Immunodeficiency Virus Type 1 RNA from Plasma samples with low viral load. Journal of Human Virology 3: 27–34, (2000)

## 7.4 Product use restriction / warranty

**NucleoSpin<sup>®</sup> RNA Virus** kit components are intended, developed, designed, and sold FOR RESEARCH PURPOSES ONLY, except, however, any other function of the product being expressly described in original MACHEREY-NAGEL product leaflets.

MACHEREY-NAGEL products are intended for GENERAL LABORATORY USE ONLY! MACHEREY-NAGEL products are suited for QUALIFIED PERSONNEL ONLY! MACHEREY-NAGEL products shall in any event only be used wearing adequate PROTECTIVE CLOTHING. For detailed information please refer to the respective Material Safety Data Sheet of the product! MACHEREY-NAGEL products shall exclusively be used in an ADEQUATE TEST ENVIRONMENT. MACHEREY-NAGEL does not assume any responsibility for damages due to improper application of our products in other fields of application. Application on the human body is STRICTLY FORBIDDEN. The respective user is liable for any and all damages resulting from such application.

DNA/RNA/PROTEIN purification products of MACHEREY-NAGEL are suitable for *IN-VITRO*-USES ONLY!

ONLY MACHEREY-NAGEL products specially labeled as IVD are also suitable for *IN-VITRO*diagnostic use. Please pay attention to the package of the product. *IN-VITRO*-diagnostic products are expressly marked as IVD on the packaging.

IF THERE IS NO IVD SIGN, THE PRODUCT SHALL NOT BE SUITABLE FOR *IN-VITRO*-DIAGNOSTIC USE!

ALL OTHER PRODUCTS NOT LABELED AS IVD ARE NOT SUITED FOR ANY CLINICAL USE (INCLUDING, BUT NOT LIMITED TO DIAGNOSTIC, THERAPEUTIC AND/OR PROGNOSTIC USE).

No claim or representations is intended for its use to identify any specific organism or for clinical use (included, but not limited to diagnostic, prognostic, therapeutic, or blood banking). It is rather in the responsibility of the user or -in any case of resale of the products -in the responsibility of the reseller to inspect and assure the use of the DNA/RNA/protein purification products of MACHEREY-NAGEL for a well-defined and specific application.

MACHEREY-NAGEL shall only be responsible for the product specifications and the performance range of MN products according to the specifications of in-house quality control, product documentation and marketing material.

This MACHEREY-NAGEL product is shipped with documentation stating specifications and other technical information. MACHEREY-NAGEL warrants to meet the stated specifications. MACHEREY-NAGEL's sole obligation and the customer's sole remedy is limited to replacement of products free of charge in the event products fail to perform as warranted. Supplementary reference is made to the general business terms and conditions of MACHEREY-NAGEL, which are printed on the price list. Please contact us if you wish to get an extra copy.

There is no warranty for and MACHEREY-NAGEL is not liable for damages or defects arising in shipping and handling (transport insurance for customers excluded), or out of accident or improper or abnormal use of this product; defects in products or components not manufactured by MACHEREY-NAGEL, or damages resulting from such non-MACHEREY-NAGEL components or products.

MACHEREY-NAGEL makes no other warranty of any kind whatsoever, and SPECIFICALLY DISCLAIMS AND EXCLUDES ALL OTHER WARRANTIES OF ANY KIND OR NATURE WHATSOEVER, DIRECTLY OR INDIRECTLY, EXPRESS OR IMPLIED, INCLUDING, WITHOUT

LIMITATION, AS TO THE SUITABILITY, REPRODUCTIVITY, DURABILITY, FITNESS FOR A PARTICULAR PURPOSE OR USE, MERCHANTABILITY, CONDITION, OR ANY OTHER MATTER WITH RESPECT TO MACHEREY-NAGEL PRODUCTS.


In no event shall MACHEREY-NAGEL be liable for claims for any other damages, whether direct, indirect, incidental, compensatory, foreseeable, consequential, or special (including but not limited to loss of use, revenue or profit), whether based upon warranty, contract, tort (including negligence) or strict liability arising in connection with the sale or the failure of MACHEREY-NAGEL products to perform in accordance with the stated specifications. This warranty is exclusive and MACHEREY-NAGEL makes no other warranty expressed or implied.

The warranty provided herein and the data, specifications and descriptions of this MACHEREY-NAGEL product appearing in MACHEREY-NAGEL published catalogues and product literature are MACHEREY-NAGEL's sole representations concerning the product and warranty. No other statements or representations, written or oral, by MACHEREY-NAGEL's employees, agent or representatives, except written statements signed by a duly authorized officer of MACHEREY-NAGEL are authorized; they should not be relied upon by the customer and are not a part of the contract of sale or of this warranty.

Product claims are subject to change. Therefore please contact our Technical Service Team for the most up-to-date information on MACHEREY-NAGEL products. You may also contact your local distributor for general scientific information. Applications mentioned in MACHEREY-NAGEL literature are provided for informational purposes only. MACHEREY-NAGEL does not warrant that all applications have been tested in MACHEREY-NAGEL laboratories using MACHEREY-NAGEL products. MACHEREY-NAGEL does not warrant the correctness of any of those applications.

Last updated: 07 / 2010, Rev. 03

Please contact: MACHEREY-NAGEL GmbH & Co. KG Tel.: +49 24 21 969-270 tech-bio@mn-net.com

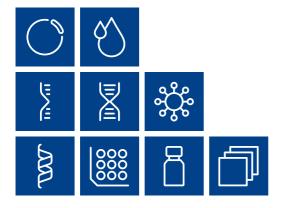

## www.mn-net.com

## **MACHEREY-NAGEL**



MACHEREY-NAGEL GmbH & Co. KG DE Tel.: +49 24 21 969-0 info@mn-net.com Valencienner Str. 11 52355 Düren · Germany

CH Tel.: +41 62 388 55 00 sales-ch@mn-net.com FR Tel.: +33 388 68 22 68 sales-fr@mn-net.com US Tel.: +1 888 321 62 24 sales-us@mn-net.com



A025752/08210