
MACHEREY-NAGEL

User manual



- Protino® GST/4B Columns 1 mL
- Protino® GST/4B Columns 5 mL

February 2022 / Rev. 01

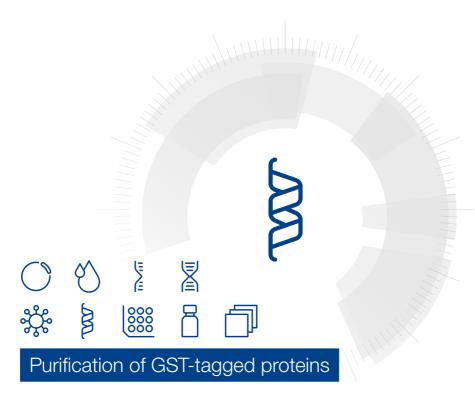

Contact MN

Germany and international

MACHEREY-NAGEL GmbH & Co. KG

Valencienner Str. 11 \cdot 52355 Düren \cdot Germany

Tel.: +49 24 21 969-0

Toll-free: 0800 26 16 000 (Germany only)

E-mail: info@mn-net.com

Technical Support Bioanalysis

Tel.: +49 24 21 969-270 E-mail: tech-bio@mn-net.com

USA

MACHEREY-NAGEL Inc.

924 Marcon Blvd. · Suite 102 · Allentown PA, 18109 · USA

Toll-free: 888 321 6224 (MACH) E-mail: sales-us@mn-net.com

France

MACHEREY-NAGEL SAS

1, rue Gutenberg - BP135 · 67720 Hoerdt Cedex · France

Tel.: +33 388 68 22 68 E-mail: sales-fr@mn-net.com

MACHEREY-NAGEL SAS (Société par Actions Simplifiée) au capital de 186600 €

Siret 379 859 531 00020 · RCS Strasbourg B379859531 · N° intracommunautaire FR04 379 859 531

Switzerland

MACHEREY-NAGEL AG

Hirsackerstr. 7 · 4702 Oensingen · Switzerland

Tel.: +41 62 388 55 00 E-mail: sales-ch@mn-net.com

www.mn-net.com

Table of contents

1	Com	ponents	4
	1.1	Kit contents and storage	4
	1.2	Reagents, consumables, and equipment to be supplied by the user	4
	1.3	About this user manual	5
2	Intro	duction	6
	2.1	Specifications	6
	2.2	Culture size	7
	2.3	Preparation of buffers	9
3	Safe	ty instructions	11
	3.1	Disposal	11
4	Prote	ocols	12
	4.1	Preparation of cleared <i>E. coli</i> lysates	12
	4.2	Purification of GST-tagged proteins using Protino® GST/4B Columns	13
5	Reg	eneration and storage	15
6	Appe	endix	16
	6.1	Troubleshooting	16
	6.2	Ordering information	18
	6.3	Product use restriction/warranty	19

1 Components

1.1 Kit contents and storage

Protino® GST/4B Columns		ıs	
	1 mL	5	mL
REF	745500.10	745515.1	745515.5
Protino® GST/4B Columns	5 x 1 mL	1 x 5 mL	5 x 5 mL
User manual	1	1	1

Shipping and storage of Protino® GST/4B Columns:

The product is shipped at ambient temperature.

Upon receipt Protino® GST/4B Columns should be **stored at 2–8 °C** (stable until: see package label). Do not freeze.

1.2 Reagents, consumables, and equipment to be supplied by the user

Reagents

- PBS (10 mM Na₂HPO₄, 1.8 mM KH₂PO₄, 2.7 mM KCl, 140 mM NaCl, pH 7.3), see section 2.3
- Elution Buffer (50 mM Tris-HCl, 10 mM glutathione, pH 8), see section 2.3
- Lysozyme (required for cell extract preparation, see section 4.1)

Consumables

Appropriate centrifuge tubes, collecting tubes

Equipment

- · Appropriate centrifuge
- Liquid chromatography system (MPLC, FPLC™, ÄKTAdesign™, etc.), peristaltic pump, or syringe
- If necessary, appropriate adaptors for connecting the Protino[®] GST/4B Columns to your system of choice. Protino[®] GST/4B Columns are equipped with 10–32 (1/16") inlet and outlet ports. With these ports the columns can easily be connected to standard MPLC/FPLC™ systems, (e.g., ÄKTAdesign™). Five adaptor sets are available for connecting the columns to other systems or for using them with a syringe (see Table 1).

Table 1: Adaptor sets					
System	Connection via	Adaptor needed	Adaptor Set		
Standard FPLC™ system (e.g., ÄKTAdesign™)	10-32	none	none		
FPLC™ system, first generation (Pharmacia)	M6	1 x M6 female to 10 – 32 male 1 x 10 – 32 female to M6 male	Protino® M6 Adaptor Set REF 745260		
MPLC system (e.g., BioLogicTM, BIO-RAD)	1/4" 28	1 x 1/4" 28 female to 10 – 32 male 1 x 10 – 32 female to 1/4" 28 female	Protino® 1/4–28 Adaptor Set REF 745261		
MPLC system (e.g., BioLogicTM, BIO-RAD)	Luer	1 x Luer female to 10 – 32 male 1 x 10 – 32 female to Luer male	Protino® Luer Adaptor Set REF 745264		
Peristaltic pump	1/16" ID tubing inlet	1 x 1/16" ID tubing to 10 – 32 male	Protino [®] Inlet PP Adaptor Set REF 745263		
Syringe	Luer inlet	1 x Luer female to 10 – 32 male	Protino [®] Inlet Luer Adaptor REF 745262		

1.3 About this user manual

It is recommended to read the instructions of this user manual carefully before use. All technical literature is also available on the internet at **www.mn-net.com**.

Please contact Technical Service regarding information about changes of the current user manual compared to previous or updated revisions.

2 Introduction

Protino® GST/4B Columns are convenient, ready-to-use FPLC™ columns prepacked with Protino® Glutathione Agarose 4B for rapid purification of GST fusion proteins. The columns can be used with an automated chromatography system, a peristaltic pump, or with a syringe for manual processing.

Protino® GST/4B Columns can be attached directly to liquid chromatography systems (such as ÄKTAdesign™ systems) via standard 10 – 32 fittings. The columns can also be operated with other chromatography systems, with a syringe, or peristaltic pump by using common adapters provided separately by MN (for details see section 1.2 or contact Technical Service Bioanalysis).

The snap-off end of the columns can be reused as stop plug for sealing the outlet of the columns for storage.

2.1 Specifications

Table 2: Specifications Protino® GST/4B Columns				
Column bed volume	1 mL	5 mL		
System compatibility	- Automated liquid chromatography systems (MPLC, FPLC™, ÄKTAdesign™, etc.) - Peristaltic pump - Syringe			
Column dimensions	0.7 cm inner diameter x 2.5 cm height	1.6 cm inner diameter x 2.5 cm height		
Column body material	Polypropylene			
Column ports	Inlet 10–32 (1/16") female Outlet 10–32 (1/16") male			
Matrix	4% beaded agarose			
Ligand	Glutathione, linked via sulfur atom			
Spacer arm	12 atoms			
Bead size	90 μm			
Binding capacity ¹	~10 mg (recombinant GST)	~50 mg (recombinant GST)		
Maximum back pressure	3 bar (0.3 MPa)			

¹ Binding capacity will vary for each GST-tagged protein.

Table 2: Specifications Protino® GST/4B Columns			
Recommended flow rates Equilibration Sample loading ² washing and elution	1.0 mL/min 2.5 mL/min 0.2 – 1.0 mL/min 0.5 – 2.0 mL/min 1.0 mL/min 5.0 mL/min		
Chemical stability	Protino® GST/4B Columns withstand incubat 0.1 M acetate pH 4, 0.1 M NaOH, 70 % ethar 6 M guanidine hydrochloride for 2 hours at ro temperature without significant loss of protein		
Storage temperature	2–8 °C 20 % ethanol		
Storage solution			

2.2 Culture size

The yield of GST-tagged proteins depends on various parameters, such as nature of the fusion protein, expression host, culture conditions, etc. However, some recommendations on protein load and culture size can be given (see Figure 1).

Culture volume requirements are based on the following assumptions:

- Protino® GST/4B Columns have a binding capacity of ~ 10 mg of recombinant GST per 1 ml bed volume.
- Typically, the expression level of GST-tagged proteins is high, ranging from 10 to 50 mg/liter of E. coli culture.
- As a starting point we recommend to use the cell lysate from a 2000–1000 mL E. coli culture per 1 mL bed volume.

² Slow flow rates are recommended for the loading step to allow maximal binding of the GST-tagged protein.

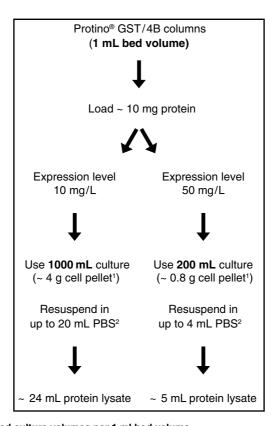


Figure 1 Required culture volumes per 1 ml bed volume

¹ On average, 250 ml of culture will produce approximately 1 g of pelleted, wet cells.

² 1 g cells may be lysed in 2-5 mL PBS, see section 4.1.

2.3 Preparation of buffers

PBS (1 liter):

10 mM	Na ₂ HPO ₄	1.780 g	Na ₂ HPO ₄ • 2 H ₂ O	$M_r = 156.01 \text{ g/mol}$	
1.8 mM	KH₂PO₄	0.245 g	KH ₂ PO ₄	$M_r = 136.09 \text{ g/mol}$	
2.7 mM	KCI	0.201 g	KCI	$M_r = 74.55 \text{ g/mol}$	
140 mM	NaCl	8.182 g	NaCl	$M_r = 58.44 \text{ g/mol}$	
Adjust pH to 7.3					

Elution Buffer (1 liter):

50 mM	Tris base	6.06 g	Tris base	$M_r = 121.14 \text{ g/mol}$
10 mM	glutathione	3.07 g	glutathione	$M_r = 307.3 \text{ g/mol}$

Adjust pH to 8.0

Prepare fresh daily and store at 4 °C

Minimum buffer volumes required for one column run can be calculated according to the following table. Note that additional volumes of PBS and Elution Buffer may be prepared to flush lines and pumps depending on the chromatographic system. As a starting point we recommend to prepare approximately 150 mL of PBS and 100 mL of Elution Buffer for the 1 mL column and 300 mL of PBS and 150 mL of Elution Buffer for the 5 mL column. Use high-purity chemicals and water for preparing the buffers. For best results, filter buffers through a 0.45 μm filter before use.

Protino® GST/4B Columns

	1 mL	5 mL
PBS		
5 mL per 1 g of cell pellet for cell extract preparation	20 mL	100 mL
10 mL per 1 mL bed volume for equilibration	10 mL	50 mL
10 mL per 1 mL bed volume for washing	10 mL	50 mL
$\sim\!50$ mL of PBS for flushing lines and pump	50 mL	50 mL
Total volume of PBS	90 mL	250 mL
Recommended volume of PBS	150 mL	300 mL

Protino® GST/4B Columns

	1 mL	5 mL
Elution Buffer		
10 mL per 1 mL bed volume for equilibration	10 mL	50 mL
$\sim\!50$ mL of Elution Buffer for flushing lines and pump	50 mL	50 mL
Total volume of Elution Buffer	60 mL	100 mL
Recommended volume of Elution Buffer	100 mL	150 mL

3 Safety instructions

When working with the **Protino**® **GST/4B Columns** kit wear suitable protective clothing (e.g., lab coat, disposable gloves, and protective goggles). For more information consult the appropriate Material Safety Data Sheets (MSDS available online at **www.mn-net.com/msds**).



The waste generated with the **Protino® GST/4B Columns** kit has not been tested for residual infectious or GMO-material. Therefore, liquid waste must be considered infectious or containing GMO-material and should be handled and discarded according local safety regulations.

3.1 Disposal

Dispose hazardous, infectious, GMO or biologically contaminated materials in a safe and acceptable manner and in accordance with all local and regulatory requirements.

4 Protocols

4.1 Preparation of cleared *E. coli* lysates

Refer to sections 2.2 for detailed information on culture volume requirements. Prepare PBS as described at section 2.3.



Refer to sections 2.2 for detailed information on culture volume requirements. Prepare PBS as described at section 2.3.

1 Cultivate and harvest cells

- As a starting point we recommend to prepare 200–1000 mL *E. coli* expression culture for the purification of 10 mg of GST-tagged protein using 1 mL bed volume using Protino[®] GST/4B Columns (see section 2.2).
- Harvest cells from an E. coli expression culture by centrifugation at 4,500–6,000 x g for 15 min at 4 °C. Remove supernatant.
- Cell pellets may be stored at -20 °C or -80 °C until needed.

2 Resuspend bacteria cells

- Thaw the cell pellet from an *E. coli* expression culture on ice (if frozen).
- Resuspend 1 g of pelleted, wet cells in 2–5 mL PBS. Pipette up and down, or stir until complete resuspension without visible cell aggregates. Perform this step on ice.

3 Lyse cells

- Add lysozyme to a final concentration of 1 mg/mL.
- · Stir the solution on ice for 30 min.
- Sonicate the suspension on ice according to the instructions provided by the manufacturer (e.g., use 10 x 15 s bursts with a 15 s cooling period between each burst).
- Carefully check samples' appearance after sonication. If the lysate is still viscous from incomplete fragmentation of DNA, add 5 μg/mL DNase I and stir on ice for 15 min.

4 Clarify lysate

- Centrifuge the crude lysate at 10,000 x g for 30 min at 4 °C to remove cellular debris.
- Carefully transfer the supernatant to a clean tube without disturbing the pellet. If the supernatant is not clear, centrifuge a second time or filter through a 0.45 μm membrane (e.g., cellulose acetate).
- Store supernatant on ice.

Proceed to section 4.2.

4.2 Purification of GST-tagged proteins using Protino® GST/4B Columns

Protino[®] GST/4B Columns can be operated with liquid chromatography systems (such as ÄKTAdesign™ systems) via standard 10–32 fittings without additional connectors.

Prepare buffers according to section 2.3. For best results, filter buffers through a 0.45 µm filter before use. Clear samples by centrifugation and/or pass them through a 0.45 µm filter.

Binding kinetics between GST and immobilized glutathione is relatively slow. Therefore use low flow rates for the loading step to allow maximal binding of the GST-tagged protein.



Prepare Elution Buffer as described at section 2.3. Elution Buffer has to be prepared fresh daily and stored at 4° C.

Protino® GST/4B columns

1 mL



1 Connect column to the chromatography system

- Purge the pump with PBS. Assure that all air is displaced.
- Remove the snap off end at the column outlet and save it for further use.
- Remove the upper plug from the column.
- Start the pump at a flow rate of approximately 0.3 mL/min.
- Fill the inlet port of the column with several drops of PBS to remove air to form a
 positive meniscus.
- Insert the fitting "drop-to-drop" into the column port to avoid introducing air bubbles.

<u>Note:</u> The snap off end can be reused as a stop plug for sealing the column outlet for storage.

2 Column equilibration

 Equilibrate the column with 5–10 column volumes of PBS until the baseline at 280 nm is stable.

5-10 mL

50-100 mL

Use a flow rate up to

1 mL/min

2.5 mL/min

3 Binding

Load the centrifuged or filtered sample onto the column.

Use a flow rate up to

0.2-1.0 mL/min

0.5-2 mL/min

<u>Note:</u> Binding kinetics between GST and immobilized glutathione is relatively slow. Therefore use low flow rates for the loading step to allow maximal binding of the GST-tagged protein.

 Collect flow through and analyze, for example by SDS-PAGE to verify that the GST-tagged protein has bound. If the fusion protein is found in early fractions of the flowthrough, the flow rate should be decreased. If the fusion protein is absent in early fractions and does appear in late fractions of the flowthrough the column capacity has been exceeded. In this case protein load should be reduced or bed volume should be increased.

4 Washing

 Wash the column with 10 column volumes of PBS or until the baseline at 280 nm is stable.

10 mL 50 mL
1 mL/min 5 mL/min

Use a flow rate up to

5 Elution

 Elute the GST-tagged protein with 10 column volumes of Elution Buffer and collect fractions.

10 mL 50 mL

1 mL/min 5 mL/min

Use a flow rate up to

 Use a Bradford protein assay, SDS-PAGE, or measure the absorbance at 280 nm to identify the fraction(s) which contain(s) the majority of the eluted GSTtagged protein and analyze by SDS-PAGE.

5 Regeneration and storage

Reuse of Protino® GST/4B Columns should only be performed with identical GST-tagged proteins to avoid possible cross-contamination. The lifetime of the matrix depends on the nature of the sample.

If a single GST-tagged protein is to be purified several times, simply wash with 10 column volumes of PBS prior to the next column run.

Basic cleaning: Wash column with approximately 10 column volumes of 100 mM Tris-HCl + 0.5 M NaCl, pH 8.5, followed by approximately 10 column volumes of 100 mM sodium acetate + 0.5 M NaCl, pH 4.5. Repeat the above wash cycles twice. Wash with 5 column volumes of PBS.

Rigorous cleaning: To remove precipitated or denatured proteins wash with 2 column volumes of 6 M guanidine hydrochloride, immediately followed by 5–10 column volumes of PBS. To remove hydrophobically bound contaminants, wash with 4 column volumes of 70 % ethanol or 1 % Triton X-100 followed by 5–10 column volumes of PBS.

If you will not be using the matrix immediately wash with additional 5 column volumes of 20 % ethanol and store at 4 $^{\circ}\text{C}$

6 Appendix

6.1 Troubleshooting

Problem

Possible cause and suggestions

Problems with vector construction

Ensure that protein and tag are in frame.

Low protein expression

Optimize bacterial expression conditions.

Low protein vield

Fusion protein forms insoluble aggregates (inclusion bodies)

Lower the growth temperature from 37 °C to 30–15 °C.

Extraction may be insufficient

- · Check extraction conditions (lysozyme, sonication).
- Use up to 2 % of a non-ionic detergent to improve cell disruption and/ or solubilization of the fusion protein.

Sonication may have been to severe

 Choose milder sonication conditions. Over-sonication can alter the conformation of the GST moiety and prevents the fusion protein from binding to Protino[®] GST/4B Columns.

Reducing agent missing

Fusion protein does not bind efficiently

 Adding DTT to the lysis buffer (final concentration 5 mM) prior to cell lysis can significantly increase binding of some fusion proteins.

Flow rate too high

 Decrease flow rate for the loading step to allow maximal binding of the GST-tagged protein.

Concentration of fusion protein is too dilute

 Concentrate the sample. Yield depends on the concentration of the fusion protein in the lysate. If the sample is too dilute, fusion proteins may not bind efficiently.

Fusion protein does not bind efficiently (continued)

Protino® GST/4B Columns have been used several times

 Clean matrix according to section 5 or use fresh matrix. Immobilized glutathione can be degraded by -glutamyl transpeptidase activity in E. coli cell lysates. Therefore, matrices with immobilized glutathione have a finite lifetime

Problem Possible cause and suggestions

Low elution volume

 Increase the volume of Elution Buffer. Depending on the nature of the fusion protein and the amount of protein loaded, additional elution steps or buffer volume is required.

Flow rate too high

Fusion protein elutes inefficiently

· Decrease flow rate during elution.

Incorrect buffer composition

 Check composition and pH of the Elution Buffer. In some cases up to 50 mM reduced glutathione may be used to improve elution.

Elution Buffer not prepared immediately before use

Prepare Elution Buffer immediately before use.

Insufficient washing

· Increase the number of washes with PBS.

Degradation of GST fusion protein

Add a protease inhibitor to the lysis solution. Multiple bands may be the result of partial degradation by host proteases during the purification procedure.

Poor protein purity

- Keep all samples and buffers on ice to reduce the activity of proteases.
- Use a protease-deficient host. Multiple bands may be the result of partial degradation by host proteases during cell growth.

Sonication may have been too severe

 Choose milder sonication conditions. Over-sonication can lead to the co-purification of host proteins with the GST-tagged protein.

Co-purification of chaperonins

Poor protein purity

(continued)

Several chaperonins, that are involved in protein folding, may copurify with GST fusion proteins, for example, DnaK (~ 70 kDa), DnaJ (~ 37 kDa), GrpE (~ 40 kDa), GroEL (~ 57 kDa), GrpE (~ 40 kDa), GroEL (57 kDa), GroES (~ 10 kDa). Several additional purification steps have been described. For example co-purification of DnaK can be avoided by treating the cell lysate with 5 mM MgCl₂ and 5 mM ATP prior to purification. DnaK can be dissociated from other components in the presence of ATP and Mg²⁺.

6.2 Ordering information

Product	REF	Pack of
Protino® GST/4B Columns 1 ml	745510.5	5 columns
Protino® GST/4B Columns 5 mL	745515.1 745515.5	1 column 5 columns
Protino® Glutathione Agarose 4B	745500.10	10 mL (settled agarose beads)
	745500.100	100 mL (settled agarose beads)
Protino® M6 Adaptor Set	745260	1 set
Protino® 1/4–28 Adaptor Set	745261	1 set
Protino® Luer Adaptor Set	745264	1 set
Protino® Inlet PP Adaptor Set	745263	1 set
Protino® Inlet Luer Adaptor	745262	1

Visit www.mn-net.com for more detailed product information.

6.3 Product use restriction/warranty

Protino® Glutathione Agarose 4B products are intended, developed, designed, and sold FOR RESEARCH PURPOSES ONLY, except, however, any other function of the product being expressly described in original MACHEREY-NAGEL product leaflets.

MACHEREY-NAGEL products are intended for GENERAL LABORATORY USE ONLY! MACHEREY-NAGEL products are suited for QUALIFIED PERSONNEL ONLY! MACHEREY-NAGEL products shall in any event only be used wearing adequate PROTECTIVE CLOTHING. For detailed information please refer to the respective Material Safety Data Sheet of the product! MACHEREY-NAGEL products shall exclusively be used in an ADEQUATE TEST ENVIRONMENT. MACHEREY-NAGEL does not assume any responsibility for damages due to improper application of our products in other fields of application. Application on the human body is STRICTLY FORBIDDEN. The respective user is liable for any and all damages resulting from such application.

DNA/RNA/PROTEIN purification products of MACHEREY-NAGEL are suitable for *IN VITRO*-USES ONLY!

ONLY MACHEREY-NAGEL products specially labeled as IVD are also suitable for *IN VITRO*-diagnostic use. Please pay attention to the package of the product. *IN VITRO*-diagnostic products are expressly marked as IVD on the packaging.

IF THERE IS NO IVD SIGN, THE PRODUCT SHALL NOT BE SUITABLE FOR *IN VITRO*-DIAGNOSTIC USE!

ALL OTHER PRODUCTS NOT LABELED AS IVD ARE NOT SUITED FOR ANY CLINICAL USE (INCLUDING, BUT NOT LIMITED TO DIAGNOSTIC, THERAPEUTIC AND/OR PROGNOSTIC USE).

No claim or representations is intended for its use to identify any specific organism or for clinical use (included, but not limited to diagnostic, prognostic, therapeutic, or blood banking). It is rather in the responsibility of the user or - in any case of resale of the products - in the responsibility of the reseller to inspect and assure the use of the DNA/RNA/protein purification products of MACHEREY-NAGEL for a well-defined and specific application.

MACHEREY-NAGEL shall only be responsible for the product specifications and the performance range of MN products according to the specifications of in-house quality control, product documentation and marketing material.

This MACHEREY-NAGEL product is shipped with documentation stating specifications and other technical information. MACHEREY-NAGEL warrants to meet the stated specifications. MACHEREY-NAGEL's sole obligation and the customer's sole remedy is limited to replacement of products free of charge in the event products fail to perform as warranted. Supplementary reference is made to the general business terms and conditions of MACHEREY-NAGEL, which are printed on the price list. Please contact us if you wish to get an extra copy.

There is no warranty for and MACHEREY-NAGEL is not liable for damages or defects arising in shipping and handling (transport insurance for customers excluded), or out of accident or improper or abnormal use of this product; defects in products or components not manufactured by MACHEREY-NAGEL, or damages resulting from such non-MACHEREY-NAGEL components or products.

MACHEREY-NAGEL makes no other warranty of any kind whatsoever, and SPECIFICALLY DISCLAIMS AND EXCLUDES ALL OTHER WARRANTIES OF ANY KIND OR NATURE WHATSOEVER, DIRECTLY OR INDIRECTLY, EXPRESS OR IMPLIED, INCLUDING,

WITHOUT LIMITATION, AS TO THE SUITABILITY, REPRODUCTIVITY, DURABILITY, FITNESS FOR A PARTICULAR PURPOSE OR USE, MERCHANTABILITY, CONDITION, OR ANY OTHER MATTER WITH RESPECT TO MACHEREY-NAGEL PRODUCTS.

In no event shall MACHEREY-NAGEL be liable for claims for any other damages, whether direct, indirect, incidental, compensatory, foreseeable, consequential, or special (including but not limited to loss of use, revenue or profit), whether based upon warranty, contract, tort (including negligence) or strict liability arising in connection with the sale or the failure of MACHEREY-NAGEL products to perform in accordance with the stated specifications. This warranty is exclusive and MACHEREY-NAGEL makes no other warranty expressed or implied.

The warranty provided herein and the data, specifications and descriptions of this MACHEREY-NAGEL product appearing in MACHEREY-NAGEL published catalogues and product literature are MACHEREY-NAGEL's sole representations concerning the product and warranty. No other statements or representations, written or oral, by MACHEREY-NAGEL's employees, agent or representatives, except written statements signed by a duly authorized officer of MACHEREY-NAGEL are authorized; they should not be relied upon by the customer and are not a part of the contract of sale or of this warranty.

Product claims are subject to change. Therefore please contact our Technical Service Team for the most up-to-date information on MACHEREY-NAGEL products. You may also contact your local distributor for general scientific information. Applications mentioned in MACHEREY-NAGEL literature are provided for informational purposes only. MACHEREY-NAGEL does not warrant that all applications have been tested in MACHEREY-NAGEL laboratories using MACHEREY-NAGEL products. MACHEREY-NAGEL does not warrant the correctness of any of those applications.

Last updated: 07 / 2010, Rev. 03

Please contact:

MACHEREY-NAGEL GmbH & Co. KG

Tel.: +49 24 21 969-270 tech-bio@mn-net.com

Trademarks:


ÄKTAdesig and FPLC are trademarks of GE Healthcare companies

BioLogic and Profinia are trademarks of BIO-RAD Laboratories, Inc.

Parafilm is a registered trademark of Bemis Company, Inc.

Protino is a registered trademark of MACHEREY-NAGEL

All used names and denotations can be brands, trademarks, or registered labels of their respective owner – also if they are not special denotation. To mention products and brands is only a kind of information (i.e., it does not offend against trademarks and brands and can not be seen as a kind of recommendation or assessment). Regarding these products or services we can not grant any guarantees regarding selection, efficiency, or operation.

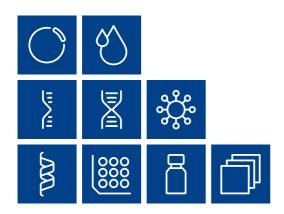

www.mn-net.com

MACHEREY-NAGEL



MACHEREY-NAGEL GmbH & Co. KG DE Tel.: +49 24 21 969-0 info@mn-net.com Valencienner Str. 11 52355 Düren · Germany

CH Tel.: +41 62 388 55 00 sales-ch@mn-net.com

FR Tel.: +33 388 68 22 68 sales-fr@mn-net.com US Tel.: +1 888 321 62 24 sales-us@mn-net.com